

## Non-destructive quantification of pharmaceutical tablet coatings using terahertz pulsed imaging and optical coherence tomography

Shuncong Zhong<sup>a,1</sup>, Yao-Chun Shen<sup>a,\*</sup>, Louise Ho<sup>b,c,2</sup>, Robert K. May<sup>d</sup>, J. Axel Zeitler<sup>d</sup>, Mike Evans<sup>c</sup>, Philip F. Taday<sup>c</sup>, Michael Pepper<sup>c</sup>, Thomas Rades<sup>b</sup>, Keith C. Gordon<sup>e</sup>, Ronny Müller<sup>f</sup>, Peter Kleinebudde<sup>f</sup>

<sup>a</sup> Department of Electrical Engineering and Electronics, University of Liverpool, L69 3GJ, UK

<sup>b</sup> School of Pharmacy, University of Otago, P.O. Box 56 Dunedin, New Zealand

<sup>c</sup> TeraView, Ltd., St. John's Innovation Park, Cambridge CB4 0WS, UK

<sup>d</sup> Department of Chemical Engineering and Biotechnology, University of Cambridge, CB2 3RA, UK

<sup>e</sup> MacDiarmid Institute for Advanced Materials and Nanotechnology, Department of Chemistry, University of Otago, P.O. Box 56 Dunedin, New Zealand

<sup>f</sup> Institute of Pharmaceutics and Biopharmaceutics, Heinrich-Heine-University, D-40225, Germany

### ARTICLE INFO

#### Article history:

Received 6 October 2010

Received in revised form

2 November 2010

Accepted 3 November 2010

Available online 3 December 2010

#### Keywords:

Non-destructive testing

Optical coherence tomography

Terahertz pulsed imaging

Tablet coating thickness

### ABSTRACT

Optical coherence tomography (OCT) and terahertz pulsed imaging (TPI) are two powerful techniques allowing high quality cross-sectional images from within scattering media to be obtained non-destructively. In this paper, we report experimental results of using OCT and TPI for quantitatively characterizing pharmaceutical tablet coatings in the thickness range of 10–140  $\mu\text{m}$ . We found that the spectral OCT system developed in-house has an axial resolution of 0.9  $\mu\text{m}$ , and is capable of quantifying very thin coatings in the range of 10–60  $\mu\text{m}$ . The upper limit of 60  $\mu\text{m}$  within the tablet coating and core is owed to the strong scattering of OCT light, which has relatively short wavelengths in the range of 0.5–1.0  $\mu\text{m}$ . On the other hand, TPI utilizes terahertz radiation that has substantially long wavelengths in the range of hundreds of microns, and thus is less prone to the scattering problem. Consequently TPI has been demonstrated to be able to quantify thicker coatings in the range of 40–140  $\mu\text{m}$  and beyond. We concluded that OCT and TPI are two complementary analytical techniques for non-destructive and quantitative characterization of pharmaceutical tablet coatings.

© 2010 Elsevier Ltd. All rights reserved.

## 1. Introduction

Pharmaceutical tablet coating is one of the preferred means to control the release of active pharmaceutical ingredients (API) molecules in the human body [1]. For example, enteric coating is used to protect the API against degradation in the stomach and sustained-release coating is used to obtain an optimum release profile, and hence a desirable API absorption rate [1]. Sustained-release coatings act as a diffusion barrier during drug release. Typically insoluble polymers are used for this purpose and depending on the polymer film thickness and porosity the drug release kinetics is controlled in such a dosage form. The quality of tablet coating will thus have direct implications on product performance. Lack of quality resulting in dose failures such as dose

dumping may precipitate legal and commercial consequences for the manufacturer. It is hence of great interest to characterize coating properties such as coating thickness and coating uniformity for the purpose of quality control and quality assurance. A number of methods have been evaluated as a means of coating quality assessment and monitoring. These methods range from coating weight gain to various destructive and non-destructive techniques [2,3]. Routinely in the pharmaceutical industry, tablet coating is controlled by employing calculations on tablet weight gain during the coating processes with respect to the amount of coating solution applied. However, weight gain determination only provides the averaged value for a batch of tablets, and it does not give information on coating thickness and its uniformity on each individual tablet [4]. Terahertz pulsed imaging (TPI) has recently been demonstrated as a powerful technique for quantitatively characterizing coatings of individual pharmaceutical tablets [5–7]. In particular, terahertz radiation can penetrate most pharmaceutical excipients, and this makes TPI an increasingly popular tool for assessing information from various pharmaceutical solid dosage forms. These include coating thickness mapping, coating interface evaluation and tablet dissolution prediction [8]. However, the thinnest coating that could be precisely quantified is about

\* Corresponding author. Tel.: +44 0 151 794 4575; fax: +44 0 151 794 4540.

E-mail address: [y.c.shen@liv.ac.uk](mailto:y.c.shen@liv.ac.uk) (Y.-C. Shen).

<sup>1</sup> Present address: School of Mechanical Engineering and Automation, Fuzhou University, Fuzhou 350108, P.R. China.

<sup>2</sup> Present address: Drug Delivery, Disposition and Dynamics, Monash Institute of Pharmaceutical Sciences, Monash University, 381 Royal Parade, Parkville, Vic. 3052, Australia.

40  $\mu\text{m}$  for most common pharmaceutical coatings, using a state-of-the-art TPI system [7,8]. Whilst this is sufficient for characterizing many finished tablet coatings, a better detection limit is highly desirable for the quality assessment of very thin tablet coatings, and for online monitoring of the coating process where such thin coatings (corresponding to the early stage of the coating process) are of particular interest.

Optical coherence tomography (OCT) is also a non-invasive and cross-sectional imaging technique that permits, for example, three-dimensional images with micrometer resolution to be obtained from within the retina [9,10]. However, most commercially available OCT systems are developed and optimized for medical applications. To date, applications of OCT beyond biomedicine have been scarce. Nevertheless, there have been some reports on using OCT for the non-destructive evaluation of thickness of papers [11], and surface roughness, gloss and bulk properties of flat-faced tablets [12]. Very recently, OCT has also been evaluated for characterizing the coating thickness of pharmaceutical tablets [13] and for imaging pharmaceutical tablets [14].

In this paper, we report an experimental investigation of using OCT and TPI for the quantitative characterization of pharmaceutical tablet coatings in the thickness range 10–140  $\mu\text{m}$ . The same set of pharmaceutical tablets was measured using both an in-house OCT system and a commercial TPI system. We found that OCT is better suited for analysis of thinner coatings owing to its high axial resolution whilst TPI is less prone to scattering problems, and thus is well positioned for characterizing pharmaceutical tablets with thicker coatings.

## 2. Methods and materials

For all TPI measurements reported here a TPI imaga 2000 system (TeraView Ltd., Cambridge, UK) was used and operated at a rapid scan mode, i.e. at an acquisition rate of 30 complete waveforms per second. At each measurement point, the terahertz radiation reflected from a tablet sample was recorded as a function of time over a scan range of a 2 mm optical delay. The TPI imaga 2000 system is specifically developed for the fully automated scan of typical pharmaceutical solid dosage forms that usually have curved surfaces. A six-axes robot system was employed to handle the tablets. This ensures that the tablet is always at the terahertz focus position with its surface perpendicular to the terahertz probe during a TPI measurement [7]. Note that the terahertz radiation used here is broadband, covering a spectral range of 5–100  $\text{cm}^{-1}$  (0.15–3 THz). The spot size of the focused terahertz beam at the tablet surface is estimated to be about 200  $\mu\text{m}$  in diameter at its centre frequency of 1.5 THz (50  $\text{cm}^{-1}$ ).

All OCT measurements were carried out using an in-house OCT system. As shown in Fig. 1, the light from a 50 W halogen lamp was

delivered onto the entrance pupil of an interferometer using a multimode fiber. The fiber end was then imaged at a ratio of 1:1 onto the surface of the sample and the reference, using a biconvex lens through a beam-splitter (50:50). The light that was reflected/scattered back from both the reference and sample was collected using another biconvex lens and then imaged back into the entrance slit of a spectrometer. Interference will occur when both the “sample image” and “reference image” are spatially matched in size and orientation within the corresponding transverse coherence distance, and their path lengths are matched within the coherence length of the light source. A PZT actuator was introduced to shift the phase of the reference beam, which allows the phase-shifting method [15] to be used to increase the measurement sensitivity. One distinct feature of our OCT system is the combination of a broadband CCD-based spectrometer with a white light source of an extremely short coherence length. This configuration provides good signal-to-noise ratio and at the same time excellent axial resolution. The multimode fiber used was chosen to have a diameter of 200  $\mu\text{m}$ , to make sure that the OCT light beam and the TPI probe beam have a comparable spot size at the sample surface.

The sample used in the present work was a batch of pharmaceutical tablets with a single coating layer. Tablet cores were biconvex (3 mm in height, 8 mm in diameter and average weight of 252 mg), and contained 10% w/w diprophyllin (API), 84.5% w/w lactose monohydrate (Flowlac<sup>®</sup> 100), 5% w/w vinylpyrrolidone-vinyl acetate copolymer (Kollidon<sup>®</sup> V64) and 0.5% w/w magnesium stearate. The coating formulation used is as follows: 50% w/w polyvinyl acetate (Kollicoat<sup>®</sup> SR 30D), 6% w/w polyvinyl alcohol-polyethyleneglycol graft copolymer (Kollicoat<sup>®</sup> IR), 0.075% w/w polyoxyethylene(20) sorbitan monooleate (Polysorbate 80), 0.3% w/w glycerolmonostearate, 0.75% w/w triethylcitrate and 42.87% w/w deionised water. The tablets were randomly selected after the following amounts of the sustained-release polymer were applied: 1.8, 3.6, 5.5, 7.3, 9.1, 10.9, 12.7, 14.5 and 18.2  $\text{mg}/\text{cm}^2$  [16].

## 3. Results and discussion

In order to characterize the performance of the in-house OCT system, we firstly measured a layered sample of the polymer film with a known thickness. Fig. 2(a) shows the broadband OCT spectral *interferogram* obtained for the layered sample (curve A), whilst curves B and C are the corresponding spectra of the sample and reference arm, respectively. The measured spectrum of the reference arm (curve C of Fig. 2(a)), which represents the overall spectral range of the instrument including the light source, has a centre wavelength of 700 nm and a full-width-at-half-maximum (FWHM) of 236 nm. The OCT axial resolution  $\Delta z$  is given by  $\Delta z = 2 \ln 2 \lambda_0^2 / \pi \Delta \lambda$ , where  $\lambda_0$  is the centre wavelength and  $\Delta \lambda$  is the FWHM of the light source. The achieved axial resolution


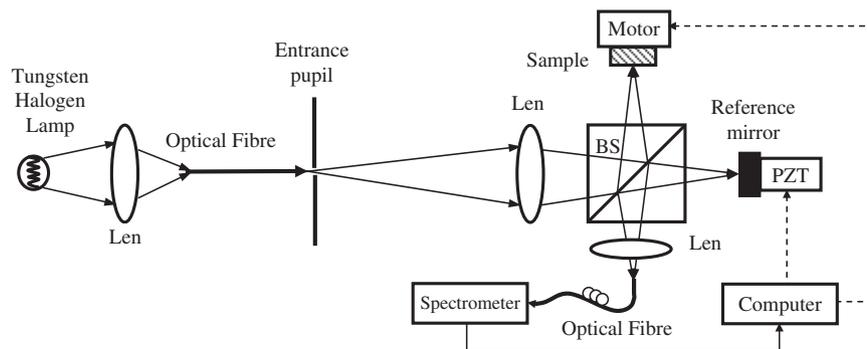


Fig. 1. Schematic diagram of spectral-domain OCT experimental setup. BS—50:50 beam-splitter; PZT—piezoelectric transducer.

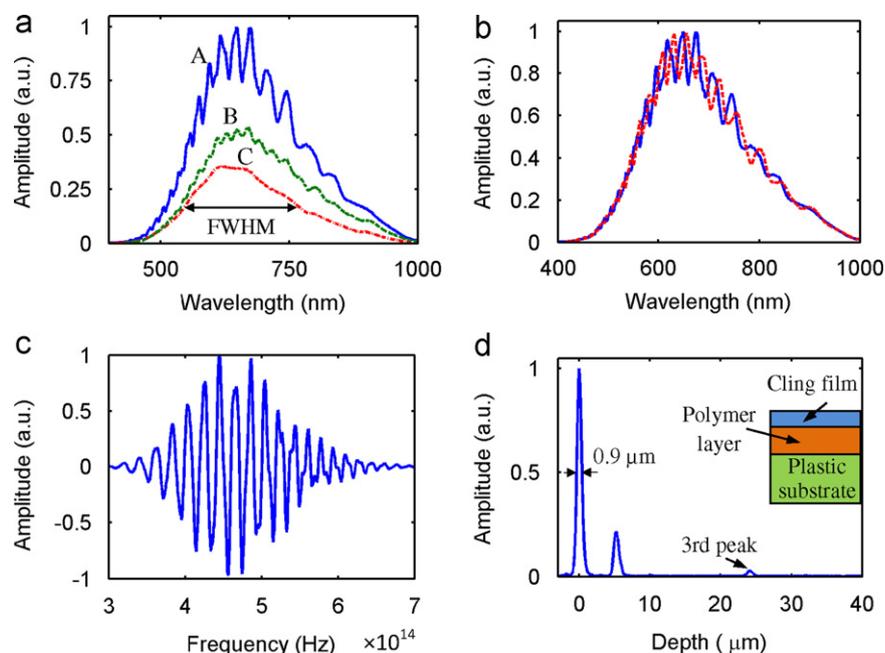

**Fig. 2.** (a) Extremely broadband OCT spectral interferogram (curve A) obtained by CCD-based interferometer. Curves B and C are the corresponding spectra measured for sample and reference arms, respectively; (b) two OCT spectral interferograms obtained with phase difference of  $\pi$ , using phase-shifting method; (c) differential OCT spectral interferogram obtained after applying phase-shift method and its corresponding OCT-FFT waveform. The axial resolution, defined as FWHM of the main peak of FFT-OCT waveform, was determined to be  $0.9 \mu\text{m}$ . Inset of (d) shows schematic diagram of a layered sample, which is a plastic plate with  $6\text{-}\mu\text{m}$ -thick cling film layer and  $18\text{-}\mu\text{m}$ -thick polymer layer.


obtained from FWHM is about  $0.9 \mu\text{m}$  in air (Fig. 2(d)), which is amongst the best ever reported [13]. Fig. 2(b) shows two spectral interferograms obtained with a phase difference of  $\pi$ , using the phase-shifting method [12]. The corresponding differential OCT interferogram (the difference of these two interferograms) is shown in Fig. 2(c). This greatly improved the signal-to-noise ratio and the measurement sensitivity. Consequently both coating layers of the sample can be resolved from its corresponding OCT-FFT waveform shown in Fig. 2(d). The thickness of the first and second layer were determined to be  $6$  and  $18 \mu\text{m}$ , respectively. Note that the third peak (at  $z=24 \mu\text{m}$ , as indicated in Fig. 2(d)) has an amplitude of  $0.027 \text{ dB}$ . This feature is very small but is still well above the background noise level ( $0.0005 \text{ dB}$ ) confirming that our OCT system has extremely good signal-to-noise ratio and thus good measurement sensitivity.

The axial resolution of the OCT setup was determined to be  $0.9 \mu\text{m}$  with a signal-to-noise ratio of 2000. A set of tablets with a film coating weight gain ranging from  $1.8$  to  $7.2 \text{ mg/cm}^2$  was measured using an OCT spectrometer. Fig. 3(a) shows a typical OCT-FFT waveform recorded for a tablet with a weight gain of  $5.4 \text{ mg/cm}^2$ . From the OCT-FFT waveform, the tablet core roughness, coating roughness and thickness can be obtained. Unlike the OCT-FFT waveform in Fig. 2(d), where a sharp reflection peak (FWHM of  $0.9 \mu\text{m}$ ) from the sample surface is observed, the OCT-FFT waveforms in Fig. 3(a) show a broad surface reflection feature comprising many sharp sparks. This can be explained as follows. (1) As discussed in the previous section the diameter of the OCT light spot size at the sample surface was chosen to be similar to that of the terahertz probe beam (about  $200 \mu\text{m}$ ) in order to make a reasonable comparison of the OCT and TPI results. (2) For a sample with a smooth surface, the surface reflection feature across the entire sampling area of  $200 \mu\text{m}$  diameter will appear at the same position ( $z=0 \mu\text{m}$ ). This leads to the sharp peak in the OCT-FFT waveform as shown in Fig. 2(c); in Fig. 3 for a rough surface, the surface height varies over the sampling spot surface. Therefore the feature reflected/scattered from different parts of the surface will


appear at a slightly different position ( $-\Delta z < z < \Delta z$ , where  $\Delta z$  is the surface roughness). This results in a broad envelope comprising many sharp sparks with each spark corresponding to light reflected/scattered from the portion of the surface of the same height.

By fitting the measured OCT-FFT waveforms with a Gaussian function, we are able to quantify the coating thickness as the peak position difference between the two fitted curves (shown as the solid curve in Fig. 3(a)). Fig. 3(b) shows the OCT b-scan maps of the same tablet with a coating polymer weight gain of  $9 \text{ mg/cm}^2$ , from which the mean thickness of this specific tablet was determined to be  $52 \mu\text{m}$ . Note that the OCT system was also used to measure tablets with a weight gain above  $9 \text{ mg/cm}^2$ . However, for these samples the OCT signal from the coating/core interface was too weak to be used for reliable coating thickness analysis. This weak OCT signal obtained from the thicker coatings was most likely caused by the strong scattering of the OCT signal within the tablet coating for this specific set of coating/core materials.

After OCT analysis, the same set of tablets was analyzed using the TPI system. Fig. 3(c) shows a typical terahertz waveform recorded for a coated tablet, with Fig. 3(d) showing the corresponding TPI b-scan map. Note here that the results shown in Fig. 3(a)–(d) were obtained using the same tablet with a weight gain of  $5.4 \text{ mg/cm}^2$ . Again, the coating thickness was calculated as the peak position difference between the surface reflection feature and the coating/core interface reflection feature. The wavelength of the terahertz radiation used for the TPI experiments is about three orders of magnitude larger than that of the OCT light. Therefore, scattering within the tablet coating, the major limiting factor in the OCT measurement, is no longer an issue in the TPI measurement. Consequently, terahertz radiation can penetrate deeper into the tablet sample and TPI is able to quantify thicker coating layers as compared to OCT. As shown in Fig. 4, TPI was successfully used for the quantitative characterization of tablet coatings of a weight gain of  $7.3\text{--}18.2 \text{ mg/cm}^2$  (thickness range of  $40\text{--}140 \mu\text{m}$ ). On the other hand, we also found that the thinnest coating layer

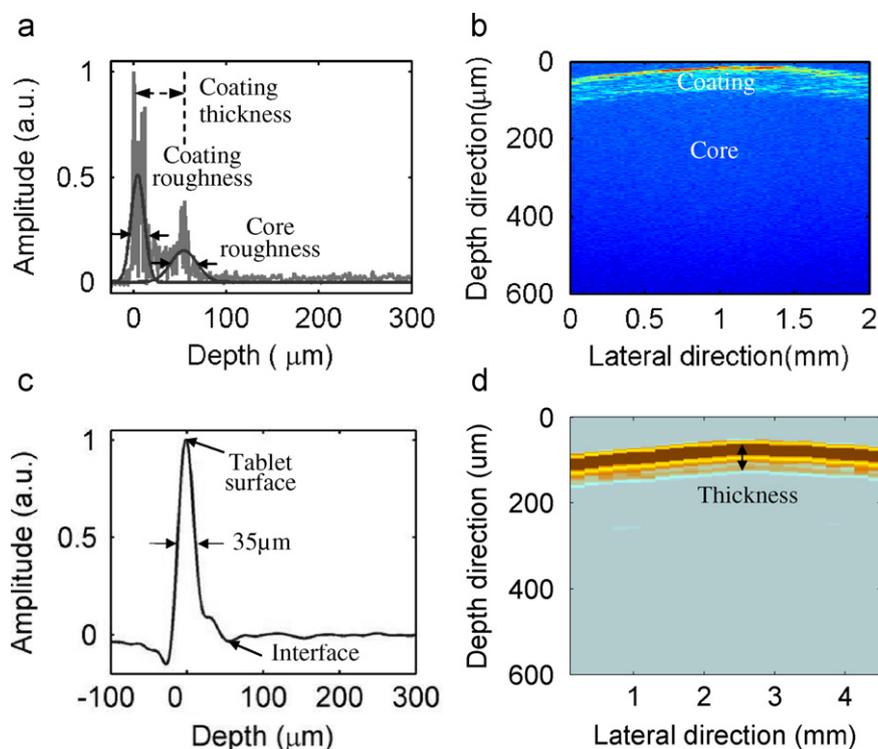

**Fig. 3.** OCT-FFT waveform (a) and OCT cross-section map (b) recorded for tablet with weight gain of  $5.4 \text{ mg/cm}^2$ . The solid lines in (a) are the best fitting curves using Gaussian function, and coating thickness was determined to be  $52 \text{ }\mu\text{m}$ . TPI terahertz waveform (c) and TPI b-scan map (d) were recorded for the same tablet. Note that TPI measures the electric field of reflected terahertz radiation and negative minimum peak at about  $52 \text{ }\mu\text{m}$  indicates that refractive index of tablet core is larger than that of tablet coating. Note here that results shown in (a)–(d) were obtained using same tablet with weight gain of  $5.4 \text{ mg/cm}^2$ .

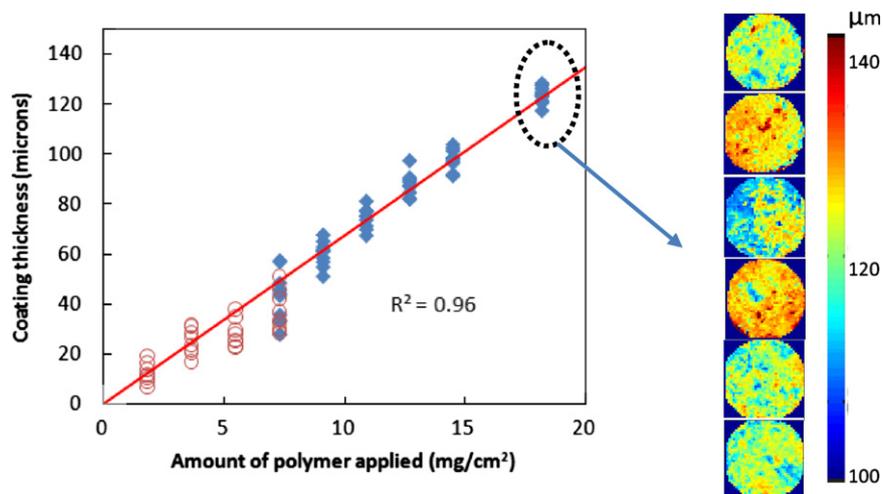

**Fig. 4.** Correlations between tablet coating thickness and amount of polymer applied. The solid diamonds represent TPI results whilst the open circles represent OCT results; inset on the right-hand side shows TPI coating thickness maps obtained for tablets with same weight gain of  $18.2 \text{ mg/cm}^2$ . There is large tablet-to-tablet variation in coating thickness for tablets with same weight gain.

that could be precisely quantified using TPI is about  $40 \text{ }\mu\text{m}$  for this specific set of tablets, which is in good agreement with previous reports [7,16]. TPI technology could not provide precise coating thickness information for tablets with weight gains of  $1.8$ ,  $3.6$  and  $5.5 \text{ mg/cm}^2$ . In contrast, OCT technology is capable of measuring such thin tablet coatings. As summarised in Fig. 4, a linear relationship ( $R^2=0.96$ ) between the coating thickness and the amount of polymer applied was observed in the full weight range of  $1.8$ – $18.2 \text{ mg/cm}^2$  (thickness range of  $10$ – $120 \text{ }\mu\text{m}$ ) investigated using both techniques. Note that neither TPI nor OCT alone can

provide coating thickness information across such a broad range of coating thickness.

#### 4. Concluding remarks

In conclusion, operating at shorter wavelength than TPI, OCT provides high axial resolution, which makes it a better choice to analyze thin coatings. Terahertz radiation has a much longer wavelength, thus being less prone to scattering problems and

hence TPI can be used to analyze thicker coatings as well as for imaging the internal structures deeper inside a tablet. The combination of OCT and TPI would provide the pharmaceutical industry with a complementary and powerful analytical technique for the quantitative characterization of pharmaceutical tablet coatings in a thickness range from a few microns to a few hundreds of microns.

### Acknowledgements

The authors would like to acknowledge the financial support from UK Technology Strategy Board (AB293H). S. Zhong would like to thank Fuzhou University Research Fund, Scientific Research Foundation for the Returned Overseas Chinese Scholars, State Education Ministry, P.R. China, and National Natural Science Foundation of China (51005077).

### References

- [1] Jantzen GM, Robinson JR. Sustained and controlled-release drug delivery systems. In: Banker GS, Rhodes CT, editors. *Modern pharmaceuticals*. 4th ed. New York: Marcel Dekker; 2002. p. 501–28.
- [2] Zeitler JA, Galdden LF. In-vitro tomography and non-destructive imaging at depth of pharmaceutical solid dosage forms. *Eur J Pharm Biopharm* 2008;71: 2–22.
- [3] Cahyadi C, Karande AD, Chan LW, Heng PWS. Comparative study of non-destructive methods to quantify thickness of tablet coatings. *Int J Pharm* 2010;398:39–49.
- [4] Perez-Ramos JD, Findlay WP, Peck G, Morris KR. Quantitative analysis of film coating in a pan coater based on in-line sensor measurements. *AAPS PharmSciTech* 2005;6(1):E127–36.
- [5] Fitzgerald AJ, Cole BE, Taday PF. Nondestructive analysis of tablet coating thickness using terahertz pulsed imaging. *J Pharm Sci* 2005;94:177–83.
- [6] Zeitler JA, Shen YC, Baker C, Taday PF, Pepper M, Rades T. Analysis of coating structures and interface in solid dosage forms by three-dimensional terahertz pulsed imaging. *J Pharm Sci* 2007;96:330–40.
- [7] Shen YC, Taday PF. Development and application of terahertz pulsed imaging for non-destructive inspection of pharmaceutical tablet. *IEEE J Sel Top Quantum Electron* 2008;14:407–15.
- [8] Ho L, Müller R, Römer M, Gordon KC, Kleinebudde P, Pepper M, Rades T, Shen YC, Taday PF, Zeitler JA. Applications of terahertz pulsed imaging to sustained-release tablet film coating quality assessment and dissolution performance. *J Control Release* 2008;127:79–87.
- [9] Maurer L, Leuenberger H. Terahertz pulsed imaging and near infrared imaging to monitor the coating process of pharmaceutical tablets. *Int J Pharm* 2009;370:8–16.
- [10] Leitgeb RA, Hitenberger CK, Fercher AF. Performance of Fourier domain vs. time domain optical coherence tomography. *Opt Express* 2003;11:889–94.
- [11] Prykari T, Czajkowski J, Alarousu E, Myllylä R. Optical coherence tomography as an accurate inspection and quality evaluation technique in paper industry. *Opt Rev* 2010;17:218–22.
- [12] Juuti M, Tuononen H, Prykari T, Kontturi V, Kuosmanen M, Alarousu E, Ketolainen J, Myllylä R, Peiponen K. Optical and terahertz measurement techniques for flat-faced pharmaceutical tablets: a case study of gloss, surface roughness and bulk properties of starch acetate tablets. *Meas Sci Technol* 2009;20:1–9.
- [13] Zhong S, Shen H, Shen YC, Zeitler JA, Ho L, Evans M, Taday PT, Pepper M, Rades T, Gordon KC, Muller R, Kleinebudde P. Noninvasive 3D characterization of layered samples using terahertz pulsed imaging and infrared optical coherence tomography. *Proceedings of the IEEE international conference on infrared, millimeter and terahertz waves*, 2009.
- [14] Mauritz JMA, Morrisby RS, Hutton RS, Legge CH, Kaminski CF. Imaging pharmaceutical tablets with optical coherence tomography. *J Pharm Sci* 2010; 99(1):385–91.
- [15] Ma ZH, Wang RK, Zhang F, Yao JQ. Spectral optical coherence tomography using two-phase shifting method. *Chin Phys Lett* 2005;22(8):1909–12.
- [16] Ho L, Müller R, Gordon KC, Kleinebudde P, Pepper M, Rades T, Shen YC, Taday PF, Zeitler JA. Terahertz pulsed imaging as an analytical tool for sustained-release tablet film coating. *Eur J Pharm Biopharm* 2009;71:117–23.